


Digital Platform KERN KDP









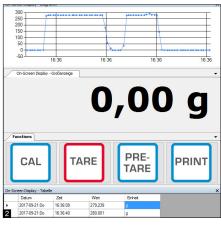

Universal weighing platforms for flexible connection in process chains

Features

- · Industry 4.0 process chains: these weighing platforms can be easily integrated into your process chain. All recorded values can be easily transferred and processed further in your system as digital data and data that has already been compared using the many interfaces available. This saves cost, time and resources and ensures even more effective working
- Connecting the platform with a tablet or PC offers the advantage of being able to use the apps and programs which are on that device. These apps and programs are usually already tailored to your needs, are easy and convenient to use and can often be expanded as needed. This offers you the maximum level of flexibility when displaying, processing and saving the recorded weighing data
- · Searching and remote control of the balance using external control devices or computers with the KERN Communication Protocol (KCP).

KCP is a standardised interface command structure for KERN balances and other instruments which allows you to recall and manage all relevant parameters and device functions. You can therefore simply connect KERN devices with KCP to computers, industrial control systems and other digital systems. In a large number of cases the KCP is compatible with the MT-SICS protocol

- · PRE-TARE function for manual subtraction of a known container weight, useful for checking fill-levels
- 1 A practical status LED informs about the active power supply
- · Draught shield standard for models with weighing plate size Ø 105 mm, weighing space W×D×H 146×146×80 mm
- · included with the delivery:
 - 2 Software BalanceConnection for adjusting and managing the KERN KDP, for large-format display of the values collected on the PC as well as transfer of this data to other Apps and programs

Technical data

- · Dimensions weighing surface
 - A Ø 105 mm, plastic, with conductive lacquer
 - B W×D 160×160 mm, stainless steel
- · Overall dimensions W×D×H 165×166×75 mm (without draught shield)
- · Cable length approx. 1,2 m
- Net weight approx. 1,2 kg
- Permissible ambient temperature 5 °C/35 °C

Accessories

- · Bluetooth data interface for wireless data transfer to PC or tablets, must be ordered at purchase, KERN KDP-A03
- WiFi data interface to transfer data from balances to PC or tablet, must be ordered at purchase, KERN YMI-A01
- Ethernet data interface, to connect an IP-based Ethernet network, continuous data transfer, must be ordered at purchase, KERN KDP-A02

Options

*WiFi and Ethernet data interface cannot be combined





















CAL EXT USB PROTOCOL	MULTI DMS 1 DAY	ET +3 DAYS	BT 4.0 WIFI LAN		
Model	Weighing capacity	Readability	Reproducibility	Linearity	Weighing plate
KERN	[Max]	[d] g	g	g	

						DAKKS Galibi. Gel tillcate
	[Max]	[d]				DAkkS
KERN	g	g	g	g		KERN
KDP 300-3	350	0,001	0,002	± 0,005	A	963-127
KDP 3000-2	3500	0,01	0,02	± 0,05	В	963-127
KDP 10K-4	10000	0,1	0,1	± 0,3	В	963-128
KDP 10K-3	10000	1	1	± 3	В	963-128



BALANCES & TEST SERVICE 2024

KERN Pictograms



Conformity Assessment

conformity assessment is

specified in the pictogram

The time required for

DAkkS calibration

DAkkS calibration

pictogram

. The time required for

is shown in days in the

The time required for

Package shipment

The time required for

in the pictogram

Pallet shipment

in the pictogram

The time required for

internal shipping prepa-

rations is shown in days

internal shipping prepa-

rations is shown in days

in days in the pictogram

Factory calibration (ISO)

Factory calibration is shown

possible (DKD)

M

DAkkS

+3 DAYS

ISO

á...



Internal adjusting

Quick setting up of the balance's accuracy with internal adjusting weight (motordriven)



Adjusting program CAL

For quick setting up of the balance's accuracy. External adjusting weight required



EasyTouch

Suitable for the connection, data transmission and control through PC or tablet



Memory

Balance memory capacity, e.g. for article data, weighing data, tare weights, PLU etc.



Alibi memory

Secure, electronic archiving of weighing results, complying with the 2014/31/EU standard.



KERN Universal Port (KUP)

allows the connection of external KUP interface adapters, e.g. RS-232, RS-485, SB, Bluetooth, WIFI, Analogue, Ethernet etc. for the exchange of data and control commands, without installation effort



RS-232 Data interface

To connect the balance to a printer, PC or network



RS-485 Data interface

To connect the balance to a printer, PC or other peripherals. Suitable for data transfer over large distances. Network in bus topology is possible



USB Data interface

To connect the balance to a printer, PC or other peripherals



Bluetooth* Data interface

To transfer data from the balance to a printer, PC or other peripherals



WIFI Data interface

To transfer data from the balance to a printer, PC or other peripherals



Control outputs

(optocoupler, digital I/O) To connect relays, signal lamps, valves, etc.



Analogue interface

to connect a suitable peripheral device for analogue processing of the measurements



Interface for second balance

For direct connection of a second balance



Network interface

For connecting the scale to an Ethernet network



KERN Communication Protocol (KCP)

It is a standardized interface command set for KERN balances and other instruments, which allows retrieving and controlling all relevant parameters and functions of the device. KERN devices featuring KCP are thus easily integrated with computers, industrial controllers and other digital systems



GLP/ISO log intern

The balance displays weight, date and time, independent of a printer connection



GLP/ISO log Printer

With weight, date and time. Only with KERN printers.



Piece counting

Reference quantities selectable. Display can be switched from piece to weight



Recipe level A

The weights of the recipe ingredients can be added together and the total weight of the recipe can be printed out



Recipe level B

Internal memory for complete recipes with name and target value of the recipe ingredients. User guidance through display



Totalising level A

The weights of similar items can be added together and the total can be printed out



Percentage determination Determining the deviation in % from the target value (100 %)



Weighing units

Can be switched to e.g. nonmetric units. See balance model. Please refer to KERN's website for more details



Weighing with tolerance range (Checkweighing)

Upper and lower limiting can be programmed individually, e.g. for sorting and dosing. The process is supported by an audible or visual signal, see the relevant model



Hold function

(Animal weighing program) When the weighing conditions are unstable, a stable weight is calculated as an average value



Protection against dust and water splashes IPxx

The type of protection is shown in the pictogram



Suspended weighing Load support with hook on the underside of the

balance



Battery operation

Ready for battery operation. The battery type is specified for each device



Rechargeable battery pack

Rechargeable set



Universal plug-in power supply

with universal input and optional input socket adapters for A) EU, CH, GB B) EU, CH, GB, US C) EU, CH, GB, US, AUS

Plug-in power supply



230V/50Hz in standard version for EU, CH. On request GB, USA or A

On request GB, USA or AUS version available



Integrated power supply unit

Integrated in balance. 230V/50Hz standard EU. More standards e.g. GB, USA or AUS on request



Weighing principle Strain gauges

Electrical resistor on an elastic deforming body



Weighing principle Tuning fork

A resonating body is electromagnetically excited, causing it to oscillate



Weighing principle Electromagnetic force compensation

Coil inside a permanent magnet. For the most accurate weighings



Weighing principle Single cell technology

Advanced version of the force compensation principle with the highest level of precision

